

ELECTROLAB

Your Quality, Our Assurance

EDI-2




Bathless Disintegration Testers


Bathless Disintegration Testers

2 Station Unit



Features:

- Bathless heating saves:
 - Time
 - Electricity 
 - Water
- Quick heating of media attains media temperature 50% faster as compared to disintegration tester with Bath
- Bathless disintegration tester enhances visual observation
- Parking of baskets at the top position at the end of the test with an alarm
- Two individual sensors measure beaker temperature
- Green illumination from bottom for better visibility of disintegration process
- Programmable timer with dual range in min-sec, hrs-min and in continuous mode
- Snap-click™ loading and a swivel free basket movement
- Maintenance free smooth drive
- Suitable for veterinary tablets and support for bolus basket

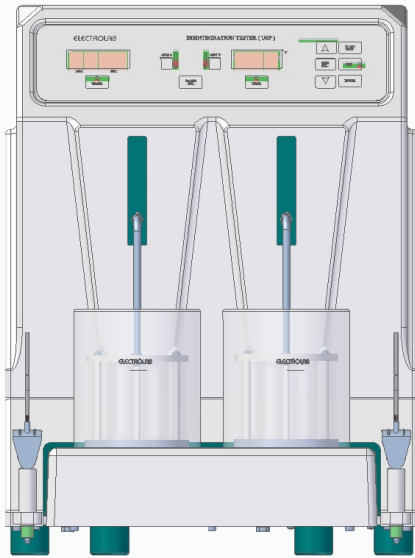

Single Time
Single Motor

Bathless Disintegration Testers

Standard Accessories



Guided disc for six tube basket assembly (Set of 6)
Part no - 0201A00004



Six tube basket assembly with guided disc (Set of 2) for EDI-2
Part no - 0202A00009



Glass beaker (1000 mL)
Part no - 0213A00001



External temperature probe for EDI-2
Part no - 0213B00001

Optional Accessories



Magnetic guided disc (Set of 12)
Part no - 0201A00007



Three tube bolus basket assembly with guided disc (Set of 2) for EDI-2
Part no - 0202A00010



Disc handling device
Part no - 0201A00018

Bathless Disintegration Testers

Specifications

| Models | EDI-2 |
|---------------------------|---|
| No. of stations | 2 |
| No. of strokes | 30 ± 1 strokes/min |
| Stroke height | 55 mm ± 2 mm |
| Temperature range | 20°C to 40°C |
| Temperature Accuracy | ± 1°C |
| Timer | Range 1: 1 sec to 99 min and 59 sec Range 2: 1 min to 9 hrs 59 min |
| Display | 4 digit LED |
| No. of motors | 2 |
| Illumination | Green Light |
| Dimensions (L x W x H) | 373 mm x 364 mm x 504 mm |
| Weight (Kg) | 18 |

Our Products

- Complete range of Dissolution Testers
- Dissolution Media Preparator
- Friability Testers
- Tablet Hardness Testers
- Electromagnetic Sieve Shakers
- Tap Density Tester
- Bulk Density Tester
- Powder Flow Tester
- Leak Testers
- Peristaltic Pumps

The information contained in this document is believed to be correct but ELECTROLAB accepts no liability for any errors and reserves the right to alter specifications without notice
November, 2017



ELECTROLAB (India) PVT. LTD. (CIN NO : U33110MH2002PTC136735)

Head Office : 401, Tirupati Udyog, I. B. Patel road, Off. Western Express highway, Goregaon (East), Mumbai - 400 063, India.

Tel : +91 - 22 - 4041 3131 • E-mail : sales@electrolabgroup.com

Factory : EL 23/24, T. T. C., Electronic Zone, M. I. D. C, Mahape, Navi Mumbai - 400 710, India.

Tel : +91 - 22 - 4161 3122 • Website : www.electrolabgroup.com

