


Flow-Through Cell Dissolution Tester Closed Loop with Auto Sampling

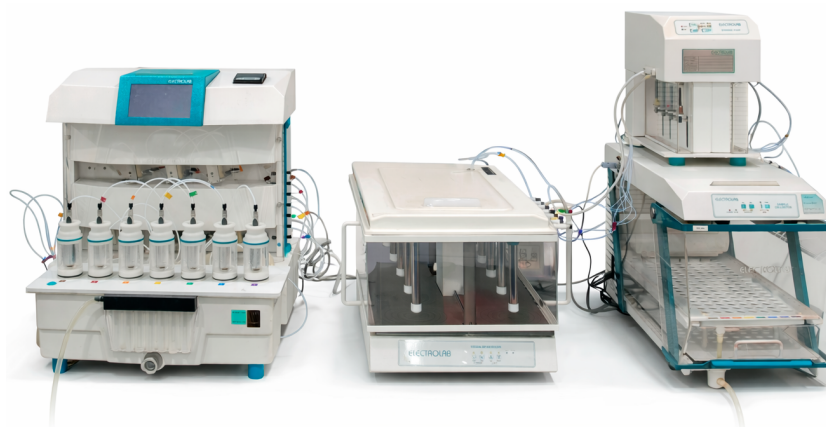

EFD-07

Features:

- Complies with USP and Ph. Eur. Specifications
 - Valveless ceramic pump heads
 - Seven flow through cell stations
 - User friendly intuitive touch screen interface
 - 1000 programmable protocols
 - Automatic flow rate adjustment for individual channels
 - Program support for gravimetric flow validation and calibration of individual channel
 - Individual cell temperature monitor
 - Isolated water circulating pump for precise
 - Temperature control of the water bath
 - Embedded printer for printing reports
- Ideal for dissolution testing of tablets, capsules, implants, suspensions, drug-eluting stents, ointments, creams, micro spheres and other sustained release dosage forms
- Visual indication of status :
Idle, Run, Error

Specifications:

Model	EFD-07 (Closed Loop with Auto Sampling)
Touch screen display	7"
Bath temperature	20°C to 45°C (Accuracy : $\pm 0.1^\circ\text{C}$) Extended temperature range from 10°C to 65°C
Temperature probe	1 no. (external)
Piston pump stroke rate	120 strokes/min (Accuracy : ± 1 stroke)
Programmable flow rate for individual channel	2 ml/min to 32 ml/min (Accuracy : $\pm 5\%$)
Dimensions (L x W x H)	1540 mm x 750 mm x 700 mm
Electrical power	220/230 VAC, 50/60 Hz, 600 W/6 Amp 110 VAC, 50/60 Hz, 600 W/6 Amp
Weight	74 Kg (Approx)

Ordering information:

Products		
Product code	Description	Quantity
1420400	EFD-07 Closed Loop with Syringe Pump and Dx Sample Collector Including Media Manager	1

Standard Accessories		
Product code	Description	Quantity
1401A00049	EFD Cell Jacket Assembly	7
1401A00058	EFD PTFE Tube Size 3.2mm OD - Length 1200mm	1
1401A00010	EFD ruby beads (5 mm)	7
1401A00002	Glass bead (1 mm)	200 gm
1401A00020	EFD Funnel	1
1401A00019	Dummy Cell Assembly (for cleaning purpose)	7

*The information contained in this document is believed to be correct but ELECTROLAB accepts no liability for any errors and reserves the right to alter specifications without notice.

ELECTROLAB (India) PVT. LTD. (CIN NO : U33110MH2002PTC136735)

Head Office: 401, Tirupati Udyog, I. B. Patel Road, Off. Western Express Highway, Goregaon (East), Mumbai - 400 063, India.

+91 - 22 - 4041 3131 sales@electrolabgroup.com

Factory: EL 23/24, T. T. C., Electronic Zone, M. I. D. C, Mahape, Navi Mumbai - 400 710, India.

+91 - 22 - 4161 3122 www.electrolabgroup.com

Flow-Through Cell Dissolution Tester

Closed Loop with Auto Sampling

Ordering information:

Optional Accessories		
Product code	Description	Quantity
1401A00001	EFD Reverse Dialysis Adapter Assembly (Set of 7) - S Type	1
1401A00008	EFD Reverse Dialysis Adapter Assembly (Set of 7)	1
1401A00009	EFD Flangeless Ferrule – Black (Set of 14)	1
1401A00014	EFD Media Selector	1
1401A00015	EFD Bath Temperature Card – 20C to 55C	1
1401A00016	EFD Flow Through Cell for Cream & Ointment (Set of 7)	1
1401A00116	EFD Media Selector Manifold	1
1401A00115	EFD Water Bath (Ver.3)	1
1401A00109	EFD Advanta Flex Tubing 4.76mm x 7.95mm For Media Selector	1
1401A00108	EFD-Scoop 22.6 mm	1
1401A00107	EFD Cell Holder	1
1401A00106	EFD Water Bath (Ver.2)	1
1401A00101	EFD-07 Front Output Manifold	1
1401A00100	EFD Close Loop Bottle Cap Assembly Including Max.Length Cannula with Flange less Nut & Ferule (Set Of 7)	1
1401A00099	EFD O Ring ID 7.2 x 1.8 Cord (Set of 7)	1
1401A00098	EFD O Ring ID 26.7 x 1.78 Cord (Set of 7)	1
1401A00097	EFD O Ring ID 6.3 x 1.8 Cord (Set of 7)	1
1401A00095	EFD O Ring ID 42.52 x 2.62 Cord (Set of 7)	1
1401A00094	EFD O Ring ID 25.12 x 1.78 Cord (Set of 7)	1
1401A00093	EFD O Ring ID 12.42 x 1.78 Cord (Set of 7)	1
1401A00092	EFD O Ring ID 26.64 x 2.62 Cord (Set of 7)	1
1401A00091	EFD O Ring ID 6.5 x 1.8 Cord (Set of 7)	1
1401A00090	EFD O Ring ID 8.92 x 1.83 Cord (Set of 7)	1
1401A000881	EFD-07-75ml Glass Vial Tray Back	1
1401A00080	EFD-07 Validation Kit	1
1401A00076	Semi Solid Adaptor for EFD (set of 7)	1
1401A00075	EFD-07 Temp Cal Cap Assy Probe	1
1401A00017	EFD Adapter for Float-A-Lyzer in 22.6 mm Cell (Set of 7)	1
1401A00071	EFD-07 Cell Kit	1
1401A00069	EFD Up Gradation From Open Loop to Close Loop	1
1401A00068	EFD-07-75ml Glass Vial Tray Front	1
1401A00067	EFD-07 Cell unit STX Comp	1
1401A00061	EFD Silicon Tubing 5 MM ID X 8 MM OD	1
1401A00060	EFD 3.2mm Ferrule (Set of 7)	1
1401A00059	EFD Flange-less nut (Set of 7)	1
1401A00018	EFD Adapter for Float-A-Lyzer in 22.6 mm Cell (Set of 7) - S Type	1
1401A00057	EFD Cannula Fixed Height for 1000ml (Set of 7)	1
1401A00056	EFD Cannula Fixed Height for 500ml (Set of 7)	1
1401A00055	EFD Cannula Fixed Height for 250ml (Set of 7)	1
1401A00054	EFD Online Connector (Set of 7)	1
1401A00053	EFD Luer Lock -BSFTL6-210-6 - Female (Set of 7)	1
1401A00052	EFD Luer Lock - BSMRLR-6 - Male (Set of 7)	1
1401A00051	EFD Luer Lock - FSLLR-1 - Female (Set of 7)	1
1401A00050	EFD Cap For Cell Jacket	1
1401A00011	EFD Quick Release Coupling With Shut of Valve (Male & Female) (Set of 2)	1
1401A00048	EFD Pressure Diverter Unit	1
1401A00013	EFD Media Manager	1
1401A00030	EFD Acrylic Spacer (Set of 7)	1

*The information contained in this document is believed to be correct but ELECTROLAB accepts no liability for any errors and reserves the right to alter specifications without notice.

ELECTROLAB (India) PVT. LTD. (CIN NO : U33110MH2002PTC136735)

Head Office: 401, Tirupati Udyog, I. B. Patel Road, Off. Western Express Highway, Goregaon (East), Mumbai - 400 063, India.

+91 - 22 - 4041 3131 sales@electrolabgroup.com

Factory: EL 23/24, T. T. C., Electronic Zone, M. I. D. C, Mahape, Navi Mumbai - 400 710, India.

+91 - 22 - 4161 3122 www.electrolabgroup.com

Flow-Through Cell Dissolution Tester

Closed Loop with Auto Sampling

Ordering information:

Optional Accessories		
Product code	Description	Quantity
1401A00029	EFD O Rings for Cell – 45mm ID (Set of 21)	1
1401A00028	EFD O Rings for Cell – 28mm ID (Set of 21)	1
1401A00024	EFD O Rings for Cell – 8.5mm ID (Set of 21)	1
1401A00023	EFD Tapper for Glass Beads (Set of 2)	1
1401A00022	EFD Sieves – 40 Mesh for 12 mm ID Powder & Granule (Set of 10)	1
1401A00021	EFD Sieves – 40 Mesh for 22.6mm Cell (Set of 10)	1

Breakdown Parts		
Product code	Description	Quantity
1401B00040	EFD-07 Platinum Kit	1
1401B00039	EFD-07 Wiring Set	1
1401B00036	EFD Cell Cap Plug Assy	1
1401B00035	EFD-Circulation Coil	1
1401B00034	EFD-07 Cell Unit Assembly	1
1401B00033	EFD-07 Solenoid Valve	1
1401B00030	EFD-07 15Pin Communication Cable (EFD to Media Selector)	1
1401B00027	EFD-07 Temperature Probe Assembly with Box	1
1401B00026	EFD-07 Individual Temperature Probe PCB with Cable	1
1401B00025	EFD-07 Piston Pump	1
1401B00024	EFD-07 Processor Card with Thermal Cut off For 55 Degree	1
1401B00023	EFD-07 Encoder Kit	1
1401B00022	EFD-07 Heater – 110V	1
1401B00021	EFD-07 Heater – 230V	1
1401B00020	EFD-07 Circulation Pump – 24V	1
1401B00019	EFD-07 Thermal Cut-off	1
1401B00018	EFD-07 Bath Temperature Probe	1
1401B00017	EFD-07 Display	1
1401B00016	EFD-07 Display Card	1
1401B00015	EFD-07 Intermediate Card	1
1401B00014	EFD-07 Printer Supply Card	1
1401B00013	EFD-07 Stepper Motor – Flow Control	1
1401B00012	EFD-07 Controller Card – Flow Rate	1
1401B00011	EFD-07 Heater Driver Card	1
1401B00010	EFD-07 SMPS – 12V	1
1401B00009	EFD-07 SMPS – 24V	1
1401B00008	EFD-07 Stepper Motor – PPM	1
1401B00007	EFD-07 Motor Opto Card – PPM	1
1401B00006	EFD-07 Motor Driver Card – PPM	1
1401B00005	EFD-07 Compatch Card	1
1401B00004	EFD-07 Software Patch	1
1401B00003	EFD-07 Processor Card with Patch	1
1401B00002	EFD-07 Driver Card for Piston Pump	1
1401B00001	EFD Individual Temperature Probe (Set of 7)	1

*The information contained in this document is believed to be correct but ELECTROLAB accepts no liability for any errors and reserves the right to alter specifications without notice.

ELECTROLAB (India) PVT. LTD. (CIN NO : U33110MH2002PTC136735)

Head Office: 401, Tirupati Udyog, I. B. Patel Road, Off. Western Express Highway, Goregaon (East), Mumbai - 400 063, India.

+91 - 22 - 4041 3131 sales@electrolabgroup.com

Factory: EL 23/24, T. T. C., Electronic Zone, M. I. D. C, Mahape, Navi Mumbai - 400 710, India.

+91 - 22 - 4161 3122 www.electrolabgroup.com

Flow-Through Cell Dissolution Tester

Closed Loop with Auto Sampling

Ordering information:

Standard Accessories for EMR		
Product code	Description	Quantity
1401A00038	EMR Vessel Kit - 1Ltr (Set of 8)	1
1401A00062	EMR Magnetic candy for 1Ltr Vessel (Set of 8)	1
1401B00032	EMR Lift Assembly (Set of 8)	1
1401A00081	EMR Tube Housing Seal (set of 8)	1
1401A00082	EMR Replenishment Vessel (1000 ml)	1
1401A00083	Replenishment Tube Housing Assembly	1

Optional Accessories for EMR		
Product code	Description	Quantity
1401A00035	EMR Vessel Kit – 60mL (Set of 7)	1
1401A00063	EMR Magnetic candy for 60mL Vessel (Set of 8)	1
1401A00085	EMR 150 ml Jar (set of 7)	1
1401A00074	150ml Bottle For EMR (Set of 7)	1
1401A00086	EMR 500 ml Vessel (Set of 7)	1
1401A00065	EMR Magnetic candy for 500mL Vessel (Set of 8)	1
1401A00066	EMR Magnetic candy for 2Ltr Vessel (Set of 8)	1
1401A00087	EMR Manual Replenishment Tube Set	1
1401A00088	EMR Magnetic candy for 150mL Vessel (Set of 8)	1
1403000289	EMR Hood Assembly	1

Assurance Plan		
Product code	Description	Quantity
1408S03	Silver Assurance Plan for Flow Through Dissolution Tester Model EFD-07 Open/Closed with Offline Accessories	1
1408S02	Gold Assurance Plan for Flow Through Dissolution Tester Model EFD-07 Open/Closed with Offline Accessories	1
1408S01	Platinum Assurance Plan for Flow Through Dissolution Tester Model EFD-07 Open/Closed with Offline Accessories	1

IQ/OQ Document		
Product code	Description	Quantity
1404S15	IQ/OQ Document for EFD-07 Closed Loop with Syringe Pump and Dx Sample Collector Including Media Manager	1

Extended IQ/OQ Document		
Product code	Description	Quantity
E1404S15	Extended IQ/OQ for EFD-07 Closed Loop with Syringe Pump and Dx Sample Collector Including Media Manager	1

Extended Warranty Support Plan		
Product code	Description	Quantity
1404S32	Extended Warranty Support Plan for EFD-07 Closed Loop with SP and Dx Sample Collector Including Media Manager	1
1404S31	Extended Warranty Support Plan for EFD-07 Open/Closed with Offline Accessories	1

Inspection & Service Charges		
Product code	Description	Quantity
1408S06	Inspection & Service Charges for EFD-07 Open/Closed with Offline Accessories	1

Re-installation Charges Without IQ/OQ		
Product code	Description	Quantity
1408S10	Re-installation Charges without IQ/OQ for EFD-07 Open/Closed with Offline Accessories	1

Single Validation Charges		
Product code	Description	Quantity
1408S04	Single Validation Charges for EFD-07 Open/Closed with Offline Accessories	1

Rental Charges		
Product code	Description	Quantity
14201R15	Rental Charges for EFD-07 Closed Loop with Manual Syringe Sampling	1

*The information contained in this document is believed to be correct but ELECTROLAB accepts no liability for any errors and reserves the right to alter specifications without notice.

ELECTROLAB (India) PVT. LTD. (CIN NO : U33110MH2002PTC136735)

Head Office: 401, Tirupati Udyog, I. B. Patel Road, Off. Western Express Highway, Goregaon (East), Mumbai - 400 063, India.

+91 - 22 - 4041 3131 sales@electrolabgroup.com

Factory: EL 23/24, T. T. C., Electronic Zone, M. I. D. C, Mahape, Navi Mumbai - 400 710, India.

+91 - 22 - 4161 3122 www.electrolabgroup.com